SUNSTAR

4200

**G·U·M**<sup>®</sup>

**PowerCARE**<sup>™</sup>  
PROFESSIONAL EXPERTISE

RECHARGEABLE  
TOOTHBRUSH

---

INTENSIVE PLAQUE  
CONTROL

---

GENTLE ON YOUR  
TEETH & GUMS

HEALTHY GUMS. HEALTHY LIFE.<sup>®</sup>

**LANGUAGES**

- 3 | **ENGLISH**
- 8 | **FRANÇAIS**
- 14 | **ITALIANO**
- 19 | **ESPAÑOL**
- 24 | **PORTUGUÊS**
- 29 | **DEUTSCH**
- 35 | **NEDERLANDS**
- 40 | **SVENSKA**
- 46 | **SUOMI**
- 51 | **POLSKI**

## READ ALL INSTRUCTIONS BEFORE USING AND SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE

### 1. INTRODUCTION

The GUM® brand builds on a basic philosophy of prevention, innovation, and quality, all on the back of SUNSTAR's pioneering vision of "sustaining a healthy body through a healthy mouth". At GUM® we are dedicated to helping people maintain healthy gums and teeth throughout their life.

#### **The Importance of Good Oral Care**

The mouth is full of bacteria. If you don't brush your teeth well or use the right products, dental plaque can build up and create tartar. Bacteria and tartar can cause an inflammation of the gums, also known as gingivitis. The gums become swollen, red, and bleed. Without early treatment, this can even develop into periodontitis which can ultimately result in tooth loss.

Power toothbrushes are more effective at plaque removal, because they help improve the brushing technique and time compliance (2 minutes).

GUM® PowerCARE™ has a unique brush head with ultra-slim bristles that provides a deep clean even in the hard to reach areas, below the gumline and between the teeth, to help maintain an optimal oral health.

### 2. PRODUCT DESCRIPTION (Figure 1)

1. Brush head with ultra-slim bristles
2. Ergonomic handle
3. On/Off/Mode button
4. Pressure sensor warning light
5. 2 Brushing modes
6. Battery light
7. Charger base with cord

### 3. SAFETY AND WARNINGS

#### **IMPORTANT SAFEGUARDS**

When using electrical products, especially when children are present, basic safety precautions should be followed, including the following:

#### **DANGER**

1. Do not immerse the charger in water or other liquids.
2. Do not use the unit while bathing.
3. Do not reach for the charger if it has fallen into water. Unplug immediately.
4. Do not place or store the charger where it can fall or be pulled into a bathtub, shower, or sink.

#### **WARNING**

To reduce the risk of burns, electrocution, fire or injury to persons:

1. Close supervision is necessary when the product is used by, on, or near children or a person with a disability.
2. This product is not intended for use by children under age three.
3. Children shall not play with the appliance.
4. Cleaning of the appliance shall not be performed by children.
5. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety

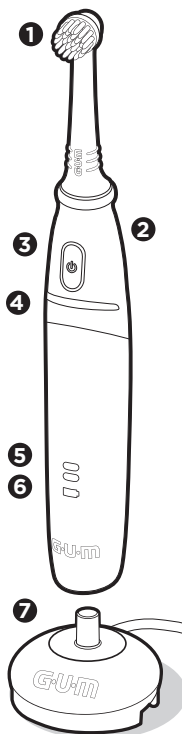
Figure 1

6. Use only for its intended use as described in this manual.
7. Never operate this product if it has a damaged cord or plug, if not working properly, if it has been dropped or damaged, or if charger has been dropped into water.
8. Keep the charger cord away from heated surfaces.
9. Never drop the appliance or insert any object into any opening of the appliance.
10. Do not use product outdoors or where aerosol (spray) products are being used. Do not use where oxygen is being administered.
11. Connect this product to a properly grounded outlet only.
12. Never force the plug into an outlet.
13. The supply cord cannot be replaced. If the cord is damaged, please discard the charger.

## 4. IMPORTANT ADDITIONAL SAFEGUARDS

Note the following for your safety:

1. Contact your dental professional before use if you have had oral or periodontal surgery in the last 12 months or have any other oral care condition currently under treatment.
2. Under normal circumstances, any light to moderate gum bleeding that occurs should stop as gums become healthy again. Contact your dental professional if excessive bleeding occurs, if bleeding continues after two weeks of brushing with this product, or if you experience any discomfort or pain when brushing.
3. Do not use a brush head with crushed or bent bristles. Replace immediately with a new brush head.
4. Periodically check the handle, brush head, charger cord, the plug, and where the cord enters the charger for any damage. Do not continue to use a damaged product.
5. Do not handle the charger or its power plug with wet hands as doing so may cause an electric shock.
6. Do not plug the charger into an electrical outlet immediately after it has been accidentally dropped into water or other liquid. Let the charger dry out for at least 24 hours before plugging into an electrical outlet.
7. If the Toothbrush is dropped, the brush head should be replaced before the next use even if no damage is visible.
8. When unplugging the charger, always do so by holding the power plug and not the cord.
9. This product contains no user-serviceable parts. Do not attempt to repair, modify or disassemble any part of the handle or charger. Doing so could cause injury, electric shock, or fire.

## 5. USING YOUR GUM® POWERCARE™ TOOTHBRUSH

### **CHARGING YOUR BRUSH**

Prior to initial use, fully charge GUM® PowerCARE™ toothbrush for at least 16 hours.

1. Plug the charger into electrical outlet. Be sure the unit is switched off before charging. It cannot be switched on during charging.
2. Full charge takes at least 16 hours and should provide 10 days of brushing based on brushing for 2 minutes twice a day.
3. When the handle is on the charger, the battery light (See Figure 1-6) will begin to flash green indicating it is charging. When the light is a steady green and no longer flashing, the handle is fully charged.
4. While using the brush, the battery light will stay green indicating the brush is fully charged.
5. When the unit needs charging, the battery light (See Figure 1-6) will turn red while the handle is on. It will then flash red for 15 seconds when the handle is turned off.

The toothbrush can be stored on the charger base (See Figure 1-7) to maintain a full charge. For the purposes of charging / recharging your GUM® PowerCARE™ toothbrush, only use the charger provided with the product.

### ASSEMBLING YOUR BRUSH

1. Insert brush head on top of the handle with bristles facing forward (see Figure 2).
2. With moderate pressure, push brush head down until bottom of the brush head stops at the top of the handle (see Figure 2).
3. If the brush head does not push down, it may not be facing forward and centered on the handle shaft. Repeat Step 1 above.

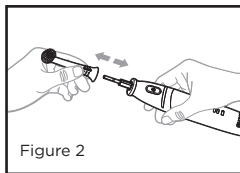
Figure 2

### PRESSURE SENSOR FOR PROTECTION FROM BRUSHING TOO HARD

1. To help protect teeth and gums from brushing too hard, a pressure sensor feature has been built into the power toothbrush.
2. If too much pressure is being applied during brushing, the pressure sensor warning light (See Figure 1-4) will turn blue and remain flashing as long as the excessive pressure continues to be applied. This is a reminder to reduce pressure while brushing. When you reduce pressure, the blue light will go off automatically.

### 2-MINUTE TIMER WITH 30-SECOND QUADRANT SIGNALS

1. Dental professionals recommend brushing the whole mouth for at least two minutes or 30 seconds for each of the four quadrants of your mouth. A quadrant includes all tooth surfaces from the front center of the mouth to the last tooth in the back. There are four quadrants, left and right, bottom and top.
2. At the end of 30, 60 and 90 seconds, there will be 2 short pulses. Those pulses signal that it is time to switch to a new quadrant. At the end of 2 minutes, there will be 4 short pulses. The brush will continue to run after 2 minutes with pulses at 30 second intervals.

### HOW TO BRUSH WITH GUM® POWERCARE™ TOOTHBRUSH

1. Apply toothpaste in a pea sized bead on the surface of the bristles.
2. Insert the toothbrush in your mouth and position the brush head on a back tooth. Press On/Off/Mode button. (See Figure 1-3) Brush in the Daily Mode or press the On/Off/Mode button again to switch to Sensitive Mode
3. For an optimal cleaning gently brush along the gumline, tooth by tooth, on both outer and inner surfaces. Then clean the chewing surface using horizontal movements.
4. Do not press too hard while brushing. At the end of each 30 -second interval, switch to another quadrant and repeat as indicated in Step 3, until the 2-minute recommended brushing time has been reached.

**Note:** Slight bleeding may occur during the first few days of use. If bleeding persists for more than two weeks, contact your dental professional.

### 2 BRUSHING MODES

1. Daily - For a deep and gentle everyday cleaning.
2. Sensitive- Extra gentle, for those who have sensitive teeth or tender gums.

### HOW TO CHANGE BRUSHING MODES

1. When the brush is turned on, (See Figure 1-3), it will automatically start in the Daily Mode. The top Mode Indicator Light (See Figure 1-5) will light up green.
2. To switch to Sensitive Mode, press On/Off /Mode button a second time. The bottom Mode Indicator Light (See Figure 1-5) will then light up green.
3. Press On/Off/Mode Button again to turn toothbrush off or twice to return to the Daily Mode.

### WHEN TO REPLACE THE BRUSH HEAD

The brush head should be replaced every 3 months as recommended by dental professionals or earlier if it is showing visible bristle wear. A worn brush head will not provide effective teeth cleaning.

### HOW TO REPLACE THE TOOTHBRUSH HEAD

1. Grip handle in one hand and brush head in the other hand.
2. Pull brush head upwards off the handle (See Figure 2).
3. Discard used brush head and insert new brush head on the handle.

4. Push the brush head down until the bottom of the brush head stops at the top of the handle. (See Figure 1-1 and 1-2).

Compatible only with GUM® PowerCARE™ replacement brush heads. Other brush heads are not approved for use on this product and may void the warranty.

### HOW TO CLEAN THE TOOTHBRUSH AND CHARGER BASE

**Important:** Clean the brush handle and brush head after every use to optimize the useful life of your power toothbrush.

#### 1. Handle and brush head:

- Rinse the brush head under running water until clean. (See Figure 3)
- With handle switched off, remove brush head. (See Figure 2) Rinse handle, including the area around the metal shaft where the brush head attaches, in lukewarm running water. (See Figure 4)
- Wipe handle and brush head dry. (See Figure 5)
- Reinsert brush head onto handle.

#### 2. Charger base:

- Unplug before cleaning.
- Clean charger base with a damp cloth. (**DO NOT** immerse in water or other liquids, rinse with water, or insert into water to clean.) Make sure that it is completely dry before plugging it back into the electrical outlet. (See Figures 6 and 7)

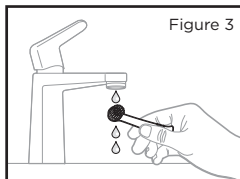
Figure 3

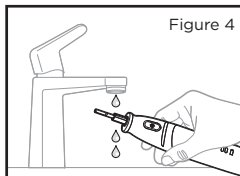
Figure 4

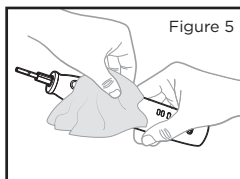
Figure 5

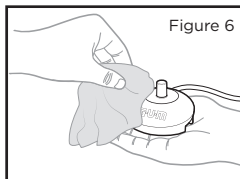
Figure 6

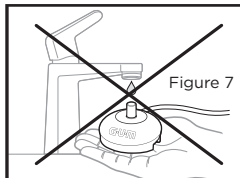
Figure 7

## 6. TECHNICAL SPECIFICATIONS

### Battery specifications

Two nickel-metal hydride (Ni-MH) rechargeable batteries: 2.4V DC

Free from cadmium, lead and mercury.

### Charger power supply

- Input: 220-240V-50Hz
- Power: 0.9W

Do not plug into an electrical outlet with a voltage other than specified above.

## 7. DISPOSAL

### Environmental Notice

This product contains recyclable electric waste and nickel-metal hydride (Ni-MH) rechargeable batteries. To protect the environment and human health, do not dispose of product in daily household waste. Once toothbrush has reached the end of its useful life, or has been damaged, please take the entire handle to your local recycling center for disposal.

Ni-MH

## 8. TROUBLESHOOTING

PROBLEM	POSSIBLE CAUSE	SOLUTION
Toothbrush does not start	<ul style="list-style-type: none"> <li>• It may require charging</li> </ul>	<ul style="list-style-type: none"> <li>• Firmly push On/Off/ Mode switch again</li> <li>• Charge the toothbrush as instructed in this manual</li> </ul>
	<ul style="list-style-type: none"> <li>• The bristles could be clogged with toothpaste.</li> </ul>	<ul style="list-style-type: none"> <li>• If the brush head bristles are clogged with toothpaste, follow the cleaning instructions in this manual</li> </ul>
Toothbrush loses its charge too quickly	<ul style="list-style-type: none"> <li>• For the first use, or after a long period of not using, the battery needs 16 hours of charge.</li> </ul>	<ul style="list-style-type: none"> <li>• Use lighter pressure when brushing.</li> <li>• Charge continuously for at least 16 hours.</li> </ul>
Brush head bristles are wearing out before the 3-month replacement period	You could be brushing too hard.	Use lighter pressure when brushing
Tickling sensation while brushing	This is normal, especially if you just switched from a manual toothbrush.	After a week of use, you will become more accustomed to the vibration.
Gums bleed slightly	<ul style="list-style-type: none"> <li>• Brushing mode is not optimal</li> <li>• May happen in the first days of usage</li> </ul>	<ul style="list-style-type: none"> <li>• Try using the "Sensitive" mode instead of the "Daily care"</li> <li>• Bleeding should stop after a few days. Should it persist more than 2 weeks, please consult your dental professional.</li> </ul>

If you cannot find a solution to your problem please contact us. You may also consult the FAQ section on [www.SunstarGUM.com](http://www.SunstarGUM.com).

## 9. WARRANTY

Sunstar Europe SA warrants this product to be free from defects in materials and workmanship for a period of two years from the date of original purchase.

Subject to the exclusions below, if the GUM® PowerCARE™ toothbrush experiences defects in materials or workmanship within two years of your original purchase date, Sunstar will replace your brush free of charge. Simply send the entire unit, including charger, postage prepaid with a copy of the original purchase receipt (price circled) and a short description of the problem.

**Exclusions:** This limited warranty does not cover brush heads, any damage caused by alterations or tampering with the product, unauthorized attempt to repair/modify the product, use contrary to the Directions for Use, normal wear and tear, purely cosmetic damage that does not affect the correct operation of the toothbrush and use of brush heads or a charger not approved for this product.

This Warranty is non-transferable. You may also have other rights which may vary country to country. This warranty is in lieu of all other warranties, express or implied or statutory, including any implied warranties of merchantability, fitness for a particular purpose or failure of essential remedy. In no event shall Sunstar Europe SA be liable for consequential or incidental damages.

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Sunstar Europe S.A., CH-1163 Etoy
Sunstar Europe S.A. Oddział w Polsce, PL-02-222 Warszawa
Sunstar France S.A.S., FR-92300 Levallois-Perret
Sunstar Deutschland GmbH., DE-79677 Schönau
Sunstar Sverige AB., SE-436 33 Askim
Sunstar Benelux B.V., NL-1382 JL Weesp