
SunstarGUM.com

(A)

(B)

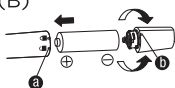4100

Made in China

4100LE126

EN DESIGNED FOR AN EXTRA CLEAN WHERE PLAQUE ACCUMULATES. PROVEN SUPERIOR EFFICACY:

Soft bristles with ultra-thin tips are 50% more efficacious in reaching the teeth and access 47x deeper below the gumline, removing plaque where a normal toothbrush can't¹. SONIC vibrations (12500/min) help remove plaque. Portable: light, cable-free, cap included. Long lasting battery - up to 90 days of use with 1 AAA battery².

Scan the QR code or visit SunstarGUM.com to learn how to best take care of your teeth and gums.

(A) HOW TO REPLACE THE BRUSH HEAD: Holding base, pull on head. Place new head on base and push down. Compatible with all GUM[®] SONIC refill heads.

(B) HOW TO REPLACE THE BATTERY: Turn and pull-down handle base. Replace it with one new AAA alkaline battery only. When replacing the battery, align symbols (a and b) on the bottom part for correct positioning. Close the container and dispose the battery according to local regulations.

HOW TO USE THE TOOTHBRUSH: Apply a small amount of toothpaste on the brush head. Insert the brush head inside the mouth. Turn on the brush pressing at the upper central button. Use the same brushing technique as with a manual toothbrush without applying any additional force.

It is recommended to change the brush head every 3 months, when bristles are worn out or after sickness.

HOW TO CLEAN THE TOOTHBRUSH: Turn off the brush. Rinse brush head with water after brushing.

WARNING! Keep reach out of children. Choking hazard. Small parts. Not for children under 3 years. For children 3 years and older, use under adult supervision only.

PACKAGING CONTENT: 1 toothbrush handle, 1 head, 1 AAA battery, 1 toothbrush cap

¹ Tested in-vitro vs. a manual flat trim toothbrush

² with 2 min brushing session twice per day.

Sunstar Europe S.A., CH-1163 Etoy
Sunstar Iberia SLU, Calle Aragón 208, 08011 Barcelona, ES

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