
SunstarGUM.com

(A)

(B)

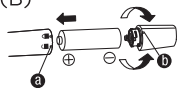DESIGNED FOR SENSITIVE TEETH AND GUMS. AN ULTRA-GENTLE & DEEP CLEAN between the teeth and under the gumline where plaque accumulates. **ULTRA-SOFT BRISTLES:** delicate on teeth and gums. **PROVEN SUPERIOR EFFICACY:** ultra-thin tip bristles remove plaque 88% deeper in-between teeth¹. SONIC vibrations (12'500/min) help remove plaque. Portable: light, cable-free, cap included. Long lasting battery – Up to 90 days of use with 1 AAA battery².

Scan the QR code or visit SunstarGUM.com to learn how to care of your sensitive teeth.

(A) HOW TO REPLACE THE BRUSH HEAD: Holding base, pull on head. Place new head on base and push down. Compatible with all GUM[®] SONIC refill heads.

(B) HOW TO REPLACE THE BATTERY: Turn and pull-down handle base. Replace it with one new AAA alkaline battery only. When replacing the battery, align symbols (a and b) on the bottom part for correct positioning. Close the container and dispose the battery according to local regulations.

HOW TO USE THE TOOTHBRUSH: Apply a small amount of toothpaste on the brush head. Insert the brush head inside the mouth. Turn on the brush pressing at the upper central button. Use the same brushing technique as with a manual toothbrush without applying any additional force.


It is recommended to change the brush head every 3 months, when bristles are worn out or after sickness.

HOW TO CLEAN THE TOOTHBRUSH: Turn off the brush. Rinse brush head with water after brushing.

WARNING! Keep reach out of children. Choking hazard. Small parts. Not for children under 3 years. For children 3 years and older, use under adult supervision only.

PACKAGING CONTENT: 1 toothbrush handle, 1 head, 1 AAA battery, 1 toothbrush cap.

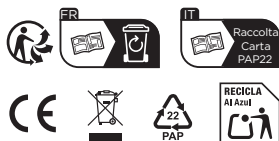¹ Tested in-vitro vs. a sensitive manual flat trim toothbrush
² with 2 min brushing session twice per day.

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Dystrybuowana przez: Sunstar Europe S.A. Oddział w Polsce, PL- 02-222, Warszawa
Distribuerad av: Sunstar Sverige AB, SE-436 33 Askim
Vertrieb durch: Sunstar Deutschland GmbH, DE-79677 Schönau
Gedistribueerd door: Sunstar Benelux B.V., NL-1382 JL Weesp